
## **Product Data Sheet**

## PE anti-mouse CD146


Catalog # / Size:	1273520 / 100 μg 1273515 / 25 μg
Clone:	ME-9F1
Isotype:	Rat IgG2a
Immunogen:	Endothelial cell line TME-3H3
<b>Reactivity:</b>	Mouse
Preparation:	The antibody was purified by affinity chromatography, and conjugated with PE under optimal conditions. The solution is free of unconjugated PE and unconjugated antibody.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
<b>Concentration:</b>	0.2

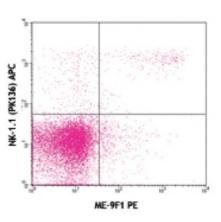

C57BL/6 splenocytes stained with ME-9F1 PE and NK-1.1 (PK136) APC

## **Applications:**

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤0.06 microg per million cells in 100 microL volume. It is recommended that the reagent be titrated for optimal performance for each application.

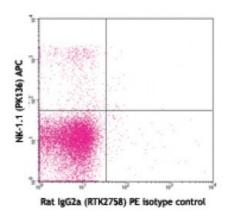

C57BL/6 splenocytes stained with rat IgG2a (RTK2758) PE isotype control and NK-1.1 (PK136) APC

Application	1. Schrage A, <i>et al.</i> 2008. <i>Histochem. Cell. Biol.</i> 129:441.
<b>References:</b>	2. Birbrair A, et al. 2013. Am J Physiol Cell Physiol. 305:1098. PubMed

**Description:** CD146, also known as melanoma cell adhesion molecule (MCAM or Mel-CAM), MUC18, S-Endo1, and A32 antigen, is an integral membrane glycoprotein that belongs to the Ig superfamily. CD146 is strongly expressed by murine vascular endothelial cells. It is expressed on about 30% of neutrophils and 60% of NK cells. Unlike in humans, CD146 is undetectable on monocytes, dendritic cells, T cells, NKT cells, B cells, or smooth muscle cells in mouse. It has been reported that an increase in CD146 expression is associated with NK cell maturation. Combined with using CD27 and CD11b staining, CD146 may be an alternative marker to detect final stages of NK cell maturation and define NK cell subsets. CD146<sup>+</sup> NK cells were found to be less cytotoxic and to produce less IFNγ than CD146- NK cells upon stimulation with target cells or activating antibodies. The role of CD146 on NK cell migration has yet to be investigated. The identification of CD146 ligand(s) will be crucial to address this issue.

Antigen 1. Despoix N, et al. 2008. Eur. J. Immunol. 38:2855.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com References: 2. Sorrentino A, *et al.* 2008. *Exp. Hematol.* 36:1035. 3. Bardin N, *et al.* 2009. *Arterioscler. Thromb. Vasc. Biol.* 29:746.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com