
Product Data Sheet

Alexa Fluor[®] 647 anti-human CD11b

| Catalog # / Size: | 2106595 / 25 tests |
|-----------------------|--|
| Clone: | ICRF44 |
| Isotype: | Mouse IgG1, κ |
| Reactivity: | Human |
| Preparation: | The antibody was purified by affinity chromatography, and conjugated with Alexa Fluor® 647 under optimal conditions. |
| Formulation: | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA). |
| Workshop Number: | IV M047 |
| Concentration: | Lot-specific |

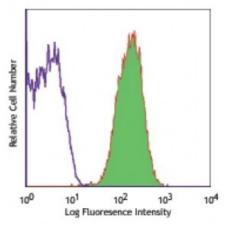

Human peripheral blood granulocytes stained with ICRF44 Alexa Fluor® 647

Applications:

Applications: Flow Cytometry

Recommended Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

* Alexa Fluor $\ensuremath{\mathbb{B}}$ 647 has a maximum emission of 668 nm when it is excited at 633nm / 635nm.

Application
Notes: The ICRF44 antibody inhibits heterotypic adhesion of granulocytes in response to fMLP. Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen tissue sections, immunofluorescence microscopy5, stimulation of monocytes3, blocking of heterotypic PMN aggregation⁸, and blocking of granulocyte activation¹². This clone was tested in-house and does not work on formalin fixed paraffin-embedded (FFPE) tissue.

The LEAF $^{\text{m}}$ purified antibody (Endotoxin <0.1 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No. 301312).

| Application | 1. Knapp W. 1989. Leucocyte Typing IV. Oxford University Press New York. |
|--------------------|---|
| References: | 2. Barclay N, <i>et al.</i> 1997. The Leucocyte Antigen Facts Book. Academic Press Inc. |
| | San Diego. |
| | 3. Rezzonico R, <i>et al.</i> 2001. <i>Blood</i> 97:2932. (Stim) |
| | 4. Marsik C, <i>et al.</i> 2003. <i>Shock</i> 20:493. (FC) |
| | 5. David A, <i>et al.</i> 2003. <i>J. Leukoc. Biol.</i> 74:551. (IF) |
| | 6. Charles N, <i>et al.</i> 2010. <i>Nat. Med.</i> 16:701. (FC) <u>PubMed</u> |
| | 7. Thurlow LR, et al. 2010. Infect. Immun. 128:1128. (FC) PubMed |
| | 8. Jadhav S, <i>et al.</i> 2001. <i>J. Immunol.</i> 167:5986. (Block) |
| | 9. Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC) |
| | 10. Sestak K, et al. 2007. Vet. Immunol. Immunopathol. 119:21. (FC) |
| | 11. Wen T, <i>et al.</i> 2014. <i>J Immunol.</i> 192:5481. (FC) <u>PubMed</u> |
| | 12. Sprong T, <i>et al.</i> 2003. <i>Blood</i> 102:3702. (Block) |
| | |

13. Cash JL, et al. 2013. EMBO Rep. 14:999. (FC) PubMed

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com **Antigen** 1. Stewart M, *et al.* 1995. *Curr. Opin. Cell Biol.* 7:690. **References:**