
## **Product Data Sheet**

## Alexa Fluor® 700 anti-human CD19

| Catalog # / Size:     | 2111130 / 100 μg<br>2111125 / 25 μg  |
|-----------------------|--|
| Clone:                | HIB19  |
| Isotype:              | Mouse IgG1, к  |
| <b>Reactivity:</b>    | Human  |
| Preparation:          | The antibody was purified by affinity chromatography, and conjugated with Alexa Fluor® 700 under optimal conditions. |
| Formulation:          | Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.  |
| Workshop<br>Number:   | V CD19.11  |
| <b>Concentration:</b> | 0.5  |

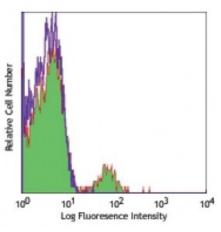

Human peripheral blood lymphocytes stained with HIB19 Alexa Fluor® 700

## **Applications:**

| Applications:              | Flow Cytometry  |
|----------------------------|---|
| Recommended<br>Usage:      | Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. The suggested use of this reagent is $\leq 1.0$ microg per million cells in 100 microL volume. It is highly recommended that the reagent be titrated for optimal performance for each application.   |
|                            | * Alexa Fluor® 700 has a maximum emission of 719 nm when it is excited at 633<br>nm / 635 nm. Prior to using Alexa Fluor® 700 conjugate for flow cytometric<br>analysis, please verify your flow cytometer's capability of exciting and detecting<br>the fluorochrome.  |
| Application<br>Notes:      | Additional reported applications (for the relevant formats) include:<br>immunohistochemical staining of acetone-fixed frozen tissue sections <sup>8</sup> and<br>blocking of B cell proliferation. Clone HIB19 is not recommended for formalin-fixed<br>paraffin-embedded sections. The LEAF <sup>™</sup> purified antibody (Endotoxin <0.1 EU/µg,<br>Azide-Free, 0.2 µm filtered) is recommended for functional assays (Cat. No.<br>302214).   |
| Application<br>References: | <ol> <li>Schlossman S, <i>et al.</i> 1995. Leucocyte Typing V. Oxford University Press. New<br/>York.</li> <li>Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV. Oxford University Press. New York.</li> <li>Bradbury L, <i>et al.</i> 1993. <i>J. Immunol.</i> 151:2915.</li> <li>Joseph A, <i>et al.</i> 2010. <i>J. Virol.</i> 84:6645. PubMed</li> <li>Wang X, <i>et al.</i> 2010. <i>Haematologica.</i> 95:884. (FC) PubMed</li> <li>Walker JD, <i>et al.</i> 2009. <i>J. Immunol.</i> 182:1548. (Block) PubMed</li> <li>Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC)</li> <li>Hansen A, <i>et al.</i> 2002. <i>Arthritis Rheum.</i> 46:2160. (IHC)</li> </ol> |
|                            |   |

**Description:** CD19 is a 95 kD type I transmembrane glycoprotein also known as B4. It is a member of the immunoglobulin superfamily expressed on B-cells (from pro-B to blastoid B cells, absent on plasma cells) and follicular dendritic cells. CD19 is involved in B cell development, activation, and differentiation. CD19 forms a complex with CD21 (CR2) and CD81 (TAPA-1), and functions as a BCR co-receptor.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com 
 Antigen
 1. Tedder T, *et al.* 1994. *Immunol. Today* 15:437.

 References:
 2. Bradbury L, *et al.* 1993. *J. Immunol.* 151:2915.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com