

**Brilliant Violet 421™ anti-human CD19**

**Catalog # / Size:** 2111170 / 100 tests  
2111165 / 25 tests

**Clone:** HIB19

**Isotype:** Mouse IgG1, κ


**Reactivity:** Human

**Preparation:** The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 421™ under optimal conditions. The solution is free of unconjugated Brilliant Violet 421™ and unconjugated antibody.

**Formulation:** Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).

**Workshop Number:** V CD19.11

**Concentration:** Lot-specific

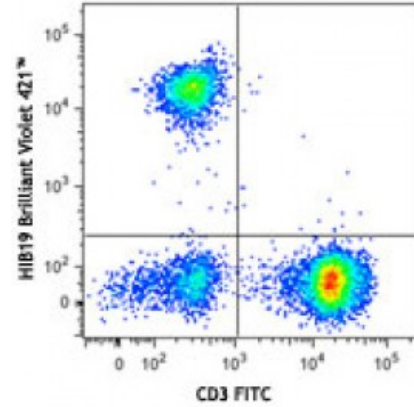


Human peripheral blood lymphocytes were stained with CD3 FITC and CD19 (clone HIB19) Brilliant Violet 421™ (top) or mouse IgG1, κ Brilliant Violet 421™ isotype control (bottom).

**Applications:**

**Applications:** Flow Cytometry

**Recommended Usage:** Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

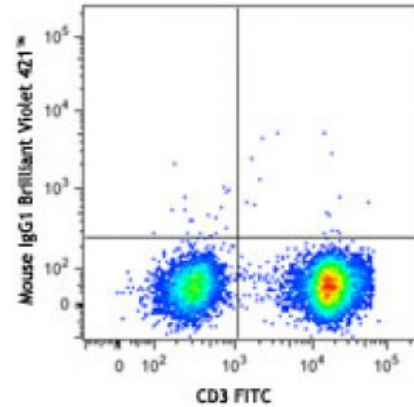

Brilliant Violet 421™ excites at 405 nm and emits at 421 nm. The standard bandpass filter 450/50 nm is recommended for detection. Brilliant Violet 421™ is a trademark of Sirigen Group Ltd.

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

**Application Notes:** Additional reported applications (for the relevant formats) include:  
immunohistochemical staining of acetone-fixed frozen tissue sections<sup>8</sup> and blocking of B cell proliferation. Clone HIB19 is not recommended for formalin-fixed paraffin-embedded sections. The LEAF™ purified antibody (Endotoxin <0.1 EU/μg, Azide-Free, 0.2 μm filtered) is recommended for functional assays (Cat. No. 302214).

**Application References:**

1. Schlossman S, *et al.* 1995. Leucocyte Typing V. Oxford University Press. New York.
2. Knapp W, *et al.* 1989. Leucocyte Typing IV. Oxford University Press. New York.
3. Bradbury L, *et al.* 1993. *J. Immunol.* 151:2915.
4. Joseph A, *et al.* 2010. *J. Virol.* 84:6645. [PubMed](#)
5. Wang X, *et al.* 2010. *Haematologica.* 95:884. (FC) [PubMed](#)
6. Walker JD, *et al.* 2009. *J. Immunol.* 182:1548. (Block) [PubMed](#)
7. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
8. Hansen A, *et al.* 2002. *Arthritis Rheum.* 46:2160. (IHC)
9. Tan YC, *et al.* 2014. *Clin Immunol.* 151:55. [PubMed](#)
10. Lu DR, *et al.* 2014. *Cell Immunol.* 152:77. [PubMed](#)

---

**Description:** CD19 is a 95 kD type I transmembrane glycoprotein also known as B4. It is a member of the immunoglobulin superfamily expressed on B-cells (from pro-B to blastoid B cells, absent on plasma cells) and follicular dendritic cells. CD19 is involved in B cell development, activation, and differentiation. CD19 forms a complex with CD21 (CR2) and CD81 (TAPA-1), and functions as a BCR co-receptor.

**Antigen References:**

1. Tedder T, *et al.* 1994. *Immunol. Today* 15:437.
2. Bradbury L, *et al.* 1993. *J. Immunol.* 151:2915.