


Brilliant Violet 421™ anti-human CD62L

Catalog # / Size:	2124140 / 100 tests 2124135 / 25 tests
Clone:	DREG-56
Isotype:	Mouse IgG1, κ
Reactivity:	Human
Preparation:	The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 421™ under optimal conditions. The solution is free of unconjugated Brilliant Violet 421™ and unconjugated antibody.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).
Workshop Number:	V S056
Concentration:	Lot-specific

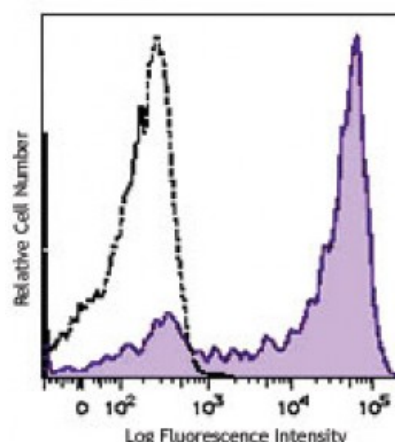

Human peripheral blood lymphocytes were stained with CD62L (clone DREG-56) Brilliant Violet 421™ (filled histogram) or mouse IgG1 Brilliant Violet 421™ isotype control (open histogram).

Applications:

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤ 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

Brilliant Violet 421™ excites at 405 nm and emits at 421 nm. The standard bandpass filter 450/50 nm is recommended for detection. Brilliant Violet 421™ is a trademark of Sirigen Group Ltd.

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Application Notes:	Additional reported applications (for the relevant formats) include: Western blotting ^{2,3,9} and <i>in vitro</i> blocking of lymphocytes binding to high endothelial venules (HEV) ² . The LEAF™ purified antibody (Endotoxin <0.1 EU/ μ g, Azide-Free, 0.2 μ m filtered) is recommended for functional assays (Cat. No. 304812).
---------------------------	--

Application References:	<ol style="list-style-type: none"> Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. Kishimoto TK, <i>et al.</i> 1990. <i>Proc. Natl. Acad. Sci. USA</i> 87:2244. (WB, Block) Jutila M, <i>et al.</i> 2002. <i>J. Immunol.</i> 169:1768. (WB) Tamassia N, <i>et al.</i> 2008. <i>J. Immunol.</i> 181:6563. (FC) PubMed Kmiecik M, <i>et al.</i> 2009. <i>J. Transl. Med.</i> 7:89. (FC) PubMed Thakral D, <i>et al.</i> 2008. <i>J. Immunol.</i> 180:7431. (FC) PubMed Charles N, <i>et al.</i> 2010. <i>Nat. Med.</i> 16:701. (FC) PubMed
--------------------------------	---

8. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
 9. Koenig JM, *et al.* 1996. *Pediatr. Res.* 39:616. (WB)
 10. Shi C, *et al.* 2011. *J. Immunol.* 187:5293. (FC) [PubMed](#)
 11. Burges M, *et al.* 2013. *Clin Cancer Res.* 19:5675. [PubMed](#)
 12. Cash JL, *et al.* 2013. *EMBO Rep.* 14:999. (FC) [PubMed](#)
-

Description: CD62L is a 74-95 kD single chain type I glycoprotein referred to as L-selectin or LECAM-1. It is expressed on most peripheral blood B cells, subsets of T and NK cells, monocytes, granulocytes, and certain hematopoietic malignant cells. CD62L binds to carbohydrates present on certain glycoforms of CD34, glycamin-1, and MAdCAM-1 and with a low affinity to anionic oligosaccharide sequences related to sialylated Lewis X (sLex, CD15s) through its C-type lectin domain. CD62L is important for the homing of naïve lymphocytes to high endothelial venules in peripheral lymph nodes and Peyer's patches. It also plays a role in leukocyte rolling on activated endothelial cells.

Antigen
References: 1. Kishimoto T, *et al.* 1990. *P. Natl. Acad. Sci. USA* 87:2244.
2. Kishimoto T, *et al.* 1991. *Blood* 78:805.