
## **Product Data Sheet**

## Alexa Fluor® 647 anti-human CD66b

Catalog # / Size:	2125545 / 25 tests 2125550 / 100 tests
Clone:	G10F5
Isotype:	Mouse IgM, κ
<b>Reactivity:</b>	Human
Preparation:	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 647 under optimal conditions.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and $0.2\%$ (w/v) BSA (origin USA).
Workshop Number:	VI MA81
Concentration:	Lot-specific

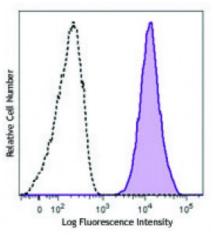

Human peripheral blood granulocytes were stained with CD66b (clone G10F5) Alexa Fluor® 647 (filled histogram) or mouse IgM, κ Alexa Fluor® 647 isotype control (open histogram).

## **Applications:**

Applications:	Flow Cytometry	and the second
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. For immunohistochemical staining on formalin-fixed paraffin-embedded tissue sections, the suggested use of this reagent is 5.0 - 10 microg per ml. It is recommended that the reagent be titrated for optimal performance for each application. * Alexa Fluor® 647 has a maximum emission of 668 nm when it is excited at 633 nm / 635 nm.	Human paraffin-embedded spleen tissue was stained with Alexa Fluor® 647 anti-human CD66b antibody (red) and co-stained with Alexa Fluor® 594 anti-human CD8a antibody (green). See additional supplemental data for detailed information.
Application Notes:	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen and formalin-fixed paraffin-embedded tissue sections.	
Application References:	<ol> <li>Schlossman S, <i>et al.</i> Eds. 1995. Leucocyt New York.</li> <li>Kishimoto T, <i>et al.</i> Eds. 1997. Leucocyte London.</li> <li>Norling LV, <i>et al.</i> 2012. Arterioscler Thro 4. Meinke P, <i>et al.</i> 2015. Neuroimmunol Dis</li> </ol>	Typing VI. Garland Publishing Inc. mb Vasc Biol. 32:1970. <u>PubMed</u>

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

lso known n pressed on to be
n oressed

 Antigen
 1. Kuijpers T, et al. 1993. J. Immunol. 151:4934.

 References:
 2. Kuroki M, et al. 1992. J. Leuk. Biol. 52:551.