
Product Data Sheet

Alexa Fluor[®] 700 anti-human TCR α/β

Catalog # / Size:	2133645 / 25 tests 2133650 / 100 tests
Clone:	IP26
Isotype:	Mouse IgG1, к
Reactivity:	Human
Preparation:	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 700 under optimal conditions.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
Concentration:	Lot-specific

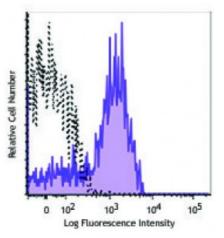

Human peripheral blood lymphocytes were stained with antihuman TCR α/β (clone IP26) Alexa Fluor® 700 (filled histogram) or mouse IgG1, κ Alexa Fluor® 700 isotype control (open histogram).

Applications:

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
	* Alexa Fluor® 700 has a maximum emission of 719 nm when it is excited at 633 nm / 635 nm. Prior to using Alexa Fluor® 700 conjugate for flow cytometric analysis, please verify your flow cytometer's capability of exciting and detecting the fluorochrome.
Application Notes:	Additional reported applications (for the relevant formats) include: T cell activation. When co-staining with anti-CD3, we recommend using clone UCHT1, since we have confirmed that IP26 does not compete with this clone. Other anti-CD3 clones may compete out the binding of IP26.
Application References:	 Schlossman S, <i>et al.</i> Eds. 1995. Leucocyte Typing V. Oxford University Press. New York. (FC) Joseph A, <i>et al.</i> 2008. <i>J. Virol.</i> 82:3078. (FC) <u>PubMed</u> Pinto JP, <i>et al.</i> 2010. <i>Immunology.</i> 130:217. <u>PubMed</u>
Description:	The IP26 antibody reacts with a monomorphic determinant of the α/β T-cell receptor, which is expressed on greater than 95% of normal peripheral blood CD3 ⁺ T cells. The α/β TCR recognizes a peptide bound to MHC leading to T-cell activation.
Antigen References:	1. Marchalonis J, <i>et al.</i> 2002. <i>J. Mol. Recognit.</i> 15:260.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com