
## **Product Data Sheet**

## Alexa Fluor® 647 anti-human Perforin

Catalog # / Size:	2140545 / 25 tests 2140550 / 100 tests
Clone:	dG9
Isotype:	Mouse lgG2b, κ
Immunogen:	Purified granules from the human lymphoma cell line
<b>Reactivity:</b>	Human
Preparation:	The antibody was purified by affinity chromatography, and conjugated with Alexa Fluor® 647 under optimal conditions.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
<b>Concentration:</b>	Lot-specific

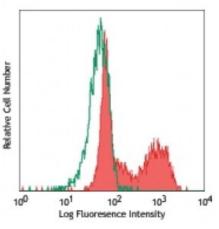

Human peripheral blood mononuclear cells intracellularly stained with DG9 Alexa Fluor® 647

## **Applications:**

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
	* Alexa Fluor ${ m I}$ 647 has a maximum emission of 668 nm when it is excited at 633nm / 635nm.
Application Notes:	Clone dG9 primarily recognizes perforin associated with cytotoxic granules <sup>9</sup> . Additional reported applications (for the relevant formats) include: immunoprecipitation, intracellular flow cytometric analysis and immunofluorescence microscopy <sup>5,7</sup> , and immunohistochemical staining of acetone-fixed frozen tissue sections and formalin-fixed paraffin-embedded tissue sections <sup>1,4</sup> .
	Does not cross-react with mouse <sup>1</sup> .
Application References:	<ol> <li>Hameed A, <i>et al.</i> 1992. <i>Am. J. Pathol.</i> 140:1025. (IHC)</li> <li>Schaerli P, <i>et al.</i> 2004. <i>J. Exp. Med.</i> 199:1265.</li> <li>Watanabe N, <i>et al.</i> 1997. <i>Blood</i> 90:3662.</li> <li>Mauad T, <i>et al.</i> 2004. <i>Pediatr. Pulmonol.</i> 38:233. (IHC)</li> <li>Barrat FJ, <i>et al.</i> 1999. <i>P. Natl. Acad. Sci. USA</i> 96:8645. (IF)</li> <li>Chen H, <i>et al.</i> 2005. <i>J. Immunol.</i> 175:591.</li> <li>Bryceson YT, <i>et al.</i> 2007. <i>Blood</i> doi:10.1182/blood-2007-02-074468. (IF)</li> <li>Wood SM, <i>et al.</i> 2009. <i>Blood</i> 114:4117. <u>PubMed</u></li> <li>Makedonas G, <i>et al.</i> 2010. <i>PLoS Pathog.</i> 6:e1000798.</li> </ol>
<b>_</b> · · ·	

**Description:** Perforin is a 70 kD cytolytic protein that is expressed in the cytoplasmic granules of cytotoxic T lymphocytes (CTLs) and natural killer (NK) cells. Perforin is one of the major effector molecules used by cytotoxic T cells and NK cells to mediate targeted cell lysis.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com Antigen1. Lieberman J. 2003. Nat. Rev. Immunol. 3:361.References:2. Trapani J, et al. 2002. Nat. Rev. Immunol. 2:735.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com