
## **Product Data Sheet**

## PerCP/Cy5.5 anti-human HLA-A,B,C

| Catalog # / Size:  | 2157100 / 100 tests<br>2157095 / 25 tests   |
|--------------------|---|
| Clone:             | W6/32   |
| Isotype:           | Mouse IgG2a, к  |
| <b>Reactivity:</b> | Human   |
| Preparation:       | The antibody was purified by affinity<br>chromatography and conjugated with<br>PerCP/Cy5.5 under optimal conditions.<br>The solution is free of unconjugated<br>PerCP/Cy5.5 and unconjugated<br>antibody. |
| Formulation:       | Phosphate-buffered solution, pH 7.2,<br>containing 0.09% sodium azide and<br>0.2% (w/v) BSA (origin USA).   |
| Concentration:     | Lot-specific  |

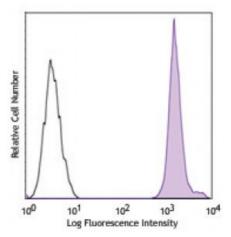

Human peripheral blood lymphocytes were stained with HLA-A,B,C (clone W6/32) PerCP/Cy5.5 (filled histogram) or mouse IgG2a PerCP/Cy5.5 isotype control (open histogram).

## **Applications:**

| Applications:              | Flow Cytometry  |
|----------------------------|---|
| Recommended<br>Usage:      | Each lot of this antibody is quality control tested by immunofluorescent staining<br>with flow cytometric analysis. For flow cytometric staining, the suggested use of<br>this reagent is 5 microL per million cells or 5 microL per 100 microL of whole<br>blood. It is recommended that the reagent be titrated for optimal performance for<br>each application.  |
|                            | * PerCP/Cy5.5 has a maximum absorption of 482 nm and a maximum emission of 690 nm.  |
| Application<br>Notes:      | Clone W6/32 recognizes a monomorphic epitope on the 45 kD polypeptide products of HLA-A, B, C <sup>18</sup> .   |
|                            | Additional reported applications (for the relevant formats) include:<br>immunoprecipitaton2, Western blotting (non-reducing)3, immunohistochemical<br>staining of acetone-fixed frozen tissue sections <sup>4,5</sup> , blocking <sup>6,7</sup> , inhibition of NK<br>cell-mediated lysis <sup>10</sup> , and activation <sup>8,9</sup> . Clone W6/32 has been reported not to be<br>suitable for immunohistochemistry on paraffin sections <sup>17</sup> . The LEAF <sup>™</sup> purified<br>antibody (Endotoxin <0.1 EU/µg, Azide-Free, 0.2 µm filtered) is recommended for<br>functional assays (Cat. No. 311412). For highly sensitive assays, we recommend<br>Ultra-LEAF <sup>™</sup> purified antibody (Cat. No. 311428) with a lower endotoxin limit than<br>standard LEAF <sup>™</sup> purified antibodies (Endotoxin <0.01 EU/microg). |
| Application<br>References: | <ol> <li>Darrow TL, <i>et al.</i> 1989. <i>J. Immunol.</i> 142:3329.</li> <li>Stern P, <i>et al.</i> 1987. <i>J. Immunol.</i> 138:1088.</li> <li>Tran TM, <i>et al.</i> 2001. <i>Immunogenetics</i> 53:440.</li> <li>Barbatis C, <i>et al.</i> 1981. <i>Gut</i> 22:985.</li> <li>Ayyoub M, <i>et al.</i> 2004. <i>Cancer Immunity</i> 4:7.</li> <li>DeFelice M, <i>et al.</i> 1990. <i>Cell. Immunol.</i> 126:420.</li> <li>Fayen J, <i>et al.</i> 1998. <i>Int. Immunol.</i> 10:1347.</li> <li>Turco MC, <i>et al.</i> 1988. <i>J. Immunol.</i> 141:2275.</li> </ol>   |

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com 9. Geppert TD, *et al.* 1989. *J. Immunol.* 142:3763.
10. Wooden SL, *et al.* 2005. *J. Immunol.* 175:1383.
11. Nagano M, *et al.* 2007. *Blood* 110:151.
12. McLoughlin RM,*et al.*2008. *J. Immunol.* 181:1323. <u>PubMed</u>
13. Takahara M, *et al.*2008. *J. Leukoc. Biol.* 83:742. <u>PubMed</u>
14. Lunemann A, *et al.*2008. *J. Immunol.* 181:6170. <u>PubMed</u>
15. Laing BJ, *et al.* 2010. *J. Thorac Cardiovasc Surg.* 139:1402. <u>PubMed</u>
16. Yoshino N, *et al.* 2000. *Exp. Anim. (Tokyo)* 49:97. (FC)
17. Vambutas A, *et al.* 2010. *J. Exp. Med.* 209:51. (epitope)

Antigen1. Barclay AN, et al. Eds. 1993. The Leukocyte Antigen FactsBook. Academic PressReferences:Inc. San Diego.