
Alexa Fluor[®] 700 anti-human TCR Vγ9

Catalog # / Size:	2256590 / 100 tests 2256585 / 25 tests	
Clone:	В3	
lsotype:	Mouse IgG1, к	001 10 ⁵
Immunogen:	PSGL-1 transfected murine 300.19 pre B-cell line	TCR Vy9 (clone 53) Alexa Fluore 700
Reactivity:	Human, Non-human primate	R Vy9 (da
Preparation:	The antibody was purified by affinity chromatography and conjugated with Alexa Fluor® 700 under optimal conditions. The solution is free of unconjugated Alexa Fluor® 700.	DT .
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).	Hur Iym CD:
Workshop Number:	V CD40.5	hur Flu к A
Concentration:	Lot-specific	con

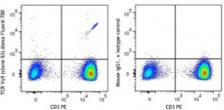

Human peripheral blood lymphocytes were stained with CD3 (clone UCHT1) PE and antihuman TCR Vγ9 (clone B3) Alexa Fluor® 700 (left) or Mouse IgG1, κ Alexa Fluor® 700 isotype control (right).

Applications:

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μ l per million cells in 100 μ l staining volume or 5 μ l per 100 μ l of whole blood.
	* Alexa Fluor® 700 has a maximum emission of 719 nm when it is excited at 633 nm / 635 nm. Prior to using Alexa Fluor® 700 conjugate for flow cytometric analysis, please verify your flow cytometer's capability of exciting and detecting the fluorochrome.
Application Notes:	Clone KPL-1 is reported to recognize the tyrosine sulfation consensus motif of PSGL-1 ¹ . Additional reported applications (for the relevant formats) include: Western Blot ¹ , immunoprecipitation ² , immunohistochemical staining of acetone-fixed frozen tissue sections and formalin-fixed paraffin embedded tissue sections ¹ , blocks the recognition of PSGL-1 with P- and L- selectin ¹ .
Application References:	1. Van Rhijn I, <i>et al.</i> 2003. <i>Intl. Immunol.</i> 15:373. 2. Yoshino N, <i>et al.</i> 2000. <i>Exp. Anim. (Tokyo)</i> 49:97. (FC)
Description:	The V γ 9 TCR is a variant of the TCR γ chain expressed on a subset of γ/δ T cells. V γ 9V δ 2 T lymphocytes, a major γ/δ T cell subset in humans, recognize phosphoantigens, certain tumor cells, and cells treated with aminobisphosphonates. This cell population displays cytolytic activity against various tumor cells. The γ/δ TCR is a heterodimeric TCR complex composed of covalently bound γ and δ chains involved in antigen recognition and the non-covalently associated monomorphic proteins CD3 δ , γ , ϵ , and ζ chains.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

Antigen	1. Scotet E, <i>et al.</i> 2005. <i>Immunity</i> 22:71
References:	2. Rincon-Orozco B, et al. 2005. J. Immunol. 175:2144

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com