

Brilliant Violet 421™ anti-human HLA-E

Catalog # / Size: 2313055 / 25 tests
2313060 / 100 tests

Clone: 3D12


Isotype: Mouse IgG1, κ

Reactivity: Human

Preparation: The antibody was purified by affinity chromatography and conjugated with Brilliant Violet 421™ under optimal conditions. The solution is free of unconjugated Brilliant Violet 421™ and unconjugated antibody.

Formulation: Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA).

Concentration: Lot-specific

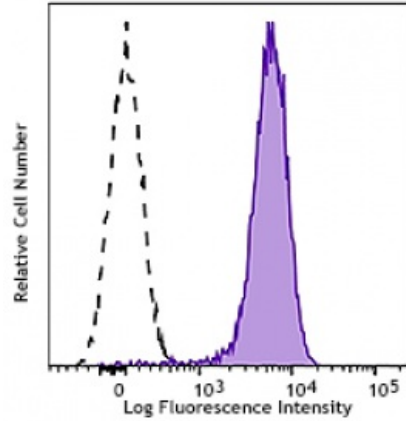

Human peripheral blood lymphocytes were stained with HLA-E (clone 3D12) Brilliant Violet 421™ (filled histogram) or mouse IgG1, κ Brilliant Violet 421™ isotype control (open histogram).

Applications:

Applications: Flow Cytometry

Recommended Usage: Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells or 5 µl per 100 µl of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.

Brilliant Violet 421™ excites at 405 nm and emits at 421 nm. The standard bandpass filter 450/50 nm is recommended for detection. Brilliant Violet 421™ is a trademark of Sirigen Group Ltd.

This product is subject to proprietary rights of Sirigen Inc. and is made and sold under license from Sirigen Inc. The purchase of this product conveys to the buyer a non-transferable right to use the purchased product for research purposes only. This product may not be resold or incorporated in any manner into another product for resale. Any use for therapeutics or diagnostics is strictly prohibited. This product is covered by U.S. Patent(s), pending patent applications and foreign equivalents.

Application References: 1. Sullivan LC, *et al.* 2008. *Tissue Antigen*. 72:415.
2. Koller BH, *et al.* 1988. *J. Immunol.* 141:897.

Description: HLA-E is a non-classical MHC class I (MHC-Ib) molecule. It is characterized by limited polymorphism and ubiquitous expression. HLA-E, as other MHC-I molecules, is heterodimerized with β2-microglobulin. HLA-E interacts with CD94/NKG2 receptor to regulate NK cell cytotoxic activities.

Antigen References: 1. Sullivan LC, *et al.* 2008. *Tissue Antigen*. 72:415.
2. Koller BH, *et al.* 1988. *J. Immunol.* 141:897.

