
Product Data Sheet

Alexa Fluor[®] 488 anti-human CD34

Catalog # / Size:	2317590 / 100 tests 2317585 / 25 tests
Clone:	581
Isotype:	Mouse IgG1, к
Reactivity:	Human
Preparation:	The antibody was purified by affinity chromatography, and conjugated with Alexa Fluor® 488 under optimal conditions.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
Workshop Number:	V MA27
Concentration:	Lot-specific

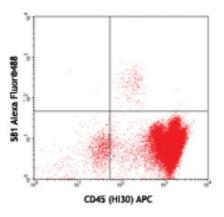

Human peripheral blood leukocytes stained with CD45 (HI30) APC and 581 Alexa Fluor® 488 (top) or mlgG1, κ Alexa Fluor® 488 isotype control (bottom). Data was generated by gating on live CD14 negative lymphocyte population.

Applications:

Applications:	Flow Cytometry	
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤ 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application. * Alexa Fluor® 488 has a maximum emission of 519 nm when it is excited at 488 nm.	
Application Notes:	The 581 antibody recognizes the class III group epitope which is resistant to sialidase/glycolyprotease and chymopapain treatment. Additional reported applications (for the relevant formats) include: immunohistochemical staining of paraffin-embedded tissue sections5 and immunofluorescence ⁶ .	
Application References:	 Schlossman SF, et al. 1995. Leukocyte Typing V:White Cell Differentiation Antigen. New York:Oxford University Press. Felschow DM, et al. 2001. Blood 97:3768. Rudin CE, et al. 1997. Br. J. Haematol. 97:488. Yoshino N, et al. 2000. Exp. Anim. (Tokyo) 49:97. (FC) Skowasch D, et al. 2003. Cardiovasc Res. 60:684. (IHC) 	

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

6. Umland O, *et al.* 2003. *J. Histochem. Cytochem.* 51:977. (IF) 7. Dong Z, *et al.* 2013. *PLoS One.* 18:80364. <u>PubMed</u>

- co co co co co co co co co co co co co c	CD34, also known as gp105-120, is a type I monomeric sialomucin-like glycophosphoprotein with an approximate molecular weight of 105-120 kD. Selectively expressed on the majority of hematopoietic stem/progenitor cells, some marrow stromal cells, capillary endothelial cells, embryonic fibroblasts, and some nervous tissue, CD34 is a commonly used marker to identify human nematopoietic stem/progenitor cells. According to the differential sensitivity to enzymatic cleavage, four groups of epitopes of CD34 have been described. CD34 mediates cell adhesion and lymphocytes homing through binding to L-selectin and E-selectin ligands.
-----------------------------------------------------------------------------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

Antigen
1. Krause DS, et al. 1996. Blood 87:1.
References:
2. Puri KD, et al. 1995. J. Cell Biol. 131:261.
3. Zola H, et al. 2007. Leukocyte and Stromal Cell Molecules: The CD Markers. John Wiley & Sons In