

Product Data Sheet

FITC anti-human CD303 (BDCA-2)

Catalog # / Size:	2371035 / 25 tests 2371040 / 100 tests
Clone:	201A
Isotype:	Mouse IgG2a, κ
Reactivity:	Human
Preparation:	The antibody was purified by affinity chromatography, and conjugated with FITC under optimal conditions. The solution is free of unconjugated FITC.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
Concentration:	Lot-specific

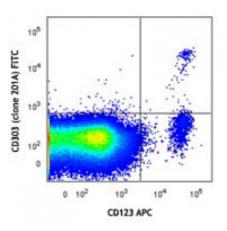

Human peripheral blood mononuclear cells were stained with CD123 APC and CD303 (clone 201A) FITC (top) or mouse IgG2a, κ FITC isotype control (bottom). Data shown was gated on the lymphocyte and monocyte populations.

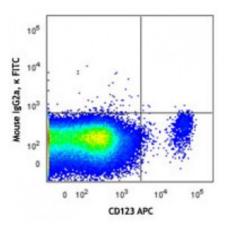

Application	1. Lee J, <i>et al.</i> 2015. <i>J Exp Med.</i> 212:385. <u>PubMed</u>
References:	2. Breton G, et al. 2015. J Exp Med. 212:401. PubMed

gl by ar	D303, also known as BDCA-2 and CLEC4C, is a 38 kD type II transmembrane lycoprotein. It is a member of the C-type lectin superfamily. CD303 is expressed y plasmacytoid dendritic cells (pDCs) and is involved in cell adhesion, signaling, nd antigen capture and processing. Crosslinking of CD303 inhibits the production f IFN- α/β and TLR-9 induced pDCs maturation.
----------------	---

Antigen	1. Jähn PS, <i>et al.</i> 2010. <i>Cell Immunol.</i> 265:15.
References:	2. Graham LM and Brown GD. 2009. <i>Cytokine</i> 48:148.
	3. Röck J, <i>et al.</i> 2007. <i>Eur. J. Immunol.</i> 37:3564.
	4. Dzionek A, <i>et al.</i> 20

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

Applications:

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.