
Product Data Sheet

PE anti-human CD21

Catalog # / Size:	2374520 / 100 tests 2374515 / 25 tests
Clone:	Bu32
Isotype:	Mouse IgG1, κ
Reactivity:	Human
Preparation:	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions. The solution is free of unconjugated PE and unconjugated antibody.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).
Workshop Number:	V CD21.4, VI CD21.5
Concentration:	Lot-specific

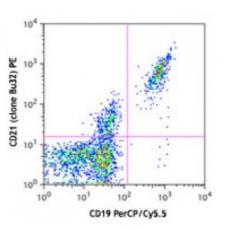

Human peripheral blood lymphocytes were stained with CD19 PerCP/Cy5.5 and CD21 (clone Bu32) PE (top) or mouse IgG1, κ PE isotype control (bottom).

104

10

Applications:

	-	<u>s</u> 1	
Applications:	Flow Cytometry	sotype contro	
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.	CD19 PerCP/Cy5.5	
Application Notes:	Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen sections4.		
Application References:	 Björck P, <i>et al.</i> 1993. <i>Eur. J. Immunol.</i> 23:1771. Frémeaux-Bacchi V, <i>et al.</i> 1996. <i>Eur. J. Immunol.</i> 26:1497. Ling NR, <i>et al.</i> 1995. <i>Clin. Exp. Immunol.</i> 101:369. Wang, C, <i>et al.</i> 2011. <i>BMC Immunol.</i> 12:53. (IHC) 		
Description:	CD21 is a 145 kD transmembrane protein also known as complement C3d receptor (C3dR), complement receptor 2 (CR2), and Epstein-Barr virus rece		

receptor (C3dR), complement receptor 2 (CR2), and Epstein-Barr virus receptor. CD21 is expressed on B cells, follicular dendritic cells, subsets of normal thymocytes and T cells, and some epithelial cells. CD21 is the receptor used by Epstein-Barr virus to infect B cells and is also the complement receptor for C3d. CD21 has also been shown to interact with a number of proteins, including CD23, CD19, annexin VI, CD81, iC3b, complement receptor 1 (CR1, CD35), and interferon- α 1 (IFN- α 1).

Antigen 1. Kishimoto T, Eds. 1997. Leukocyte Typing VI. Garland Publishing Inc.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

 Moore MD, *et al.* 1987. *Proc. Natl. Acad. Sci. USA* 84:9194.
 Szakonyi G, *et al.* 2001. *Science* 292:1725. **References:**

- 4. Weis JJ, et