
APC/Fire[™] 750 anti-human CD27

Catalog # / Size:	2382135 / 25 tests 2382140 / 100 tests	
Clone:	M-T271	Elative Cell Number
lsotype:	Mouse IgG1, к	
Immunogen:	Human T cells from a T-ALL patient.	
Reactivity:	Human, Non-human primate, Other	
Preparation:	The antibody was purified by affinity chromatography and conjugated with APC/Fire™	2 ushbby
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).	Log P
Workshop Number:	750 under optimal conditions.	Human periphe lymphocytes w CD27 (clone M
Concentration:	Lot-specific	750 (filled hist IgG1, κ APC/Fi

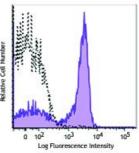

Human peripheral blood lymphocytes were stained with CD27 (clone M-T271) APC/Fire™ 750 (filled histogram) or mouse IgG1, κ APC/Fire™ 750 isotype control (open histogram).

Applications:

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 μ l per million cells in 100 μ l staining volume or 5 μ l per 100 μ l of whole blood.
	* APC/Fire™ 750 has a maximum excitation of 650 nm and a maximum emission of 787 nm.
Application Notes:	Additional reported applications (for the relevant formats) include: immunohistochemical staining of formalin-fixed paraffin-embedded frozen tissue sections ¹ , immunofluorescent staining ² , and ELISA ³ .
Application References:	 Ma S, et al. 2011. J. Virol. 85:165. (IHC) Manzo A, et al. 2008. Arthritis Rheum. 11:3377. (IF) Kato K, et al. 2007. Exp. Hematol. 35:434. (ELISA)
Description:	CD27 is a 50-55 kD type I membrane protein also known as S152 and T14. It is a lymphocyte-specific member of the TNF-receptor superfamily. CD27 is expressed on medullary thymocytes, virtually all mature T cells, some B cells, and NK cells. CD27 binds to CD70, and plays a role in costimulation of T cell activation and regulation of B cell differentiation and proliferation. The cytoplasmic domains of CD27 have also been shown to interact with TRAF2 and TRAF5 to elicit NF-kB and SAPK/JNK activation.
Antigen References:	 Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV: White Cell Differentiation Antigens. Oxford University Press. Schlossman S, <i>et al.</i> 1995. Leucocyte Typing V: White Cell Differentiation Antigens. Oxford University Press. Hintzen R, <i>et al.</i> 1994. <i>Immunol. Today</i> 15:307. Agematsu K, <i>et al.</i> 1995. <i>J. Immunol.</i> 154:3627.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com