
## PE anti-human HLA-DR, DP, DQ

Catalog # / Size:		
Clone:	Tü39	
lsotype:	Mouse IgG2a, к	H 10 <sup>6</sup>
Immunogen:	Human PBL	14 (come Tub) 19 (come Tub) 10 (co
<b>Reactivity:</b>	Human	
Preparation:	The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions. The solution is free of unconjugated PE and unconjugated antibody.	
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and 0.2% (w/v) BSA (origin USA).	Human Iymphoo
Concentration:	Lot-specific	CD19 Al DR, DP,

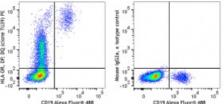

Human peripheral blood lymphocytes were stained with CD19 Alexa Fluor® 488 and HLA-DR, DP, DQ (clone TÃf¼39) PE (left) or Mouse IgG2a, κ PE isotype control (right).

## **Applications:**

Applications:	Flow Cytometry
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 $\mu$ l per million cells or 5 $\mu$ l per 100 $\mu$ l of whole blood. It is recommended that the reagent be titrated for optimal performance for each application.
Application Notes:	Tü39 has been reported to react with a shared epitope of HLA-DR, HLA-DP, and HLA-DQ.
	Additional reported applications (of relevant formats) include immunoprecipitation <sup>6</sup> , <i>in vitro</i> blocking of MLR <sup>5</sup> , and suppressor cell generation <sup>4</sup> .
Application References:	<ol> <li>Thorsby E. 1974. <i>Transplant. Rev.</i> 18:51.</li> <li>Qvigstad E, <i>et al.</i> 1984. <i>Hum. Immunol.</i> 11:207.</li> <li>Servenius B, <i>et al.</i> 1984. <i>EMBO J.</i> 3:3209.</li> <li>Ottenhoff TH, <i>et al.</i> 1985. <i>Hu</i></li> </ol>
Description:	HLA-DR, HLA-DP, and HLA-DQ are heterodimeric cell surface glycoproteins comprised of an $\alpha$ (heavy) chain and a $\beta$ (light) chain. They are expressed on B cells, activated T cells, monocytes/macrophages, dendritic cells, and

B cells, activated T cells, monocytes/macrophages, dendritic cells, and other non-professional APCs. In conjunction with the CD3/TCR complex and CD4 molecules, HLA-DR is critical for efficient peptide presentation to CD4+ T cells. Variations in the HLA gene expression are crucial to graft survival.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

Antigen	1. Thorsby E. 1974. <i>Transplant. Rev.</i> 18:51.
References:	2. Qvigstad E, et al. 1984. Hum. Immunol. 11:207.
	3. Servenius B, et al. 1984. EMBO J. 3:3209.

- Ottenhoff TH, et al. 1985. Hum. Immunol. 13:105.
   Strassmann G, et al. 1985. Hum. Immunol. 13:125.
- 6. Trowsdale J, et al. 1985. Immunol Rev. 85:5.