Product Data Sheet

APC/Fire™ 750 anti-human CD11a/CD18 (LFA-1)

Catalog # / 2417100 / 100 tests

Size: 2417095 / 25 tests

Clone: m24

Isotype: Mouse IgG1, κ

Immunogen: Fibronectin-purified human

monocytes

Reactivity: Human

Preparation: The antibody was purified by affinity

chromatography and conjugated with

APC/Fire™ 750 under optimal

conditions.


Formulation: Phosphate-buffered solution, pH 7.2,

containing 0.09% sodium azide and

0.2% (w/v) BSA (origin USA).

Workshop Number: 750 under optimal conditions.

Concentration: Lot-specific

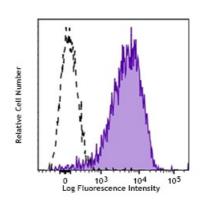

PMA-stimulated human peripheral blood granulocytes were stained with CD11a/CD18 (clone m24) APC/Fire™ 750 (filled histogram) or mouse IgG1, κ APC/Fire™ 750 isotype control (open histogram).

Applications:

Applications: Flow Cytometry

Recommended

Usage:

Each lot of this antibody is quality control tested by immunofluorescent

staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 µl per million cells in 100 µl staining volume or 5 µl per

 $100~\mu l$ of whole blood.

* APC/Fire™ 750 has a maximum excitation of 650 nm and a maximum

emission of 787 nm.

Application Notes:

Clone m24 can be used as a reporter of the activation state of integrin

receptor in response to exposure

with Mg^{2+} or Mn^{2+} .

Application References:

1. Hogg N. et al. 1985. Cell Immunol. 92:247.

2. McDowall A. et al. 2003. J. Clin. Invest. 111:51.

Description:

CD11/CD18 belongs to the integrin family of proteins. It is heterodimeric cell surface receptor expressed on all leukocytes. CD18, in association with integrin α chain CD11a, CD11b, and CD11c forms LFA-1, Mac-1, and $\alpha_X\beta_2$, respectively, and plays an important role in leukocyte adhesion. CD11/CD18 complexes bind ICAM-1 (CD54), ICAM-2 (CD102), ICAM-3 (CD50), iC3b, and fibrinogen. Clone m24 binds the extended/open high affinity conformation of CD11a/CD18. The antibody can be used as a reporter of the activation state of the integrin receptor in response to exposure to Mg²+ or Mn²+.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products.

Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

 Anderson D, et al. 1987. Annu. Rev. Med. 38:175.
 Springer T. 1994. Cell 76:301. References:

Antigen