Product Data Sheet

Alexa Fluor® 647 anti-BrdU

Catalog # / Size: 2420540 / 100 tests

2420535 / 25 tests

Clone: 3D4

Isotype: Mouse IgG1, κ

Immunogen: Iodouridine-conjugated ovalbumin

Preparation: The antibody was purified by affinity

chromatography and conjugated with

Alexa Fluor® 647 under optimal


conditions.

Formulation: Phosphate-buffered solution, pH 7.2,

containing 0.09% sodium azide and

0.2% (w/v) BSA (origin USA).

Concentration: Lot-specific

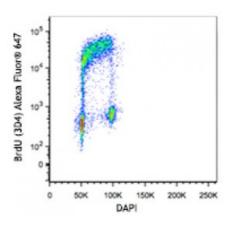

HEL cell line was pulsed with BrdU for 1 hour (upper panel) or without (lower panel) and then stained with anti-BrdU (3D4) Alexa Fluor® 647 according to BioLegend BrdU staining procedure. Cells were subsequently stained with 1 microg of DAPI for DNA

150K

DAPI

200K

Applications:

Applications: Flow Cytometry

Recommended

Usage:

Each lot of this antibody is quality control tested by intracellular immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is 5 microL per million cells or 5 microL per 100 microL of whole blood. It is recommended that the reagent be titrated for optimal performance for

each application.

* Alexa Fluor\$ 647 has a maximum emission of 668 nm when it is excited at

633 nm / 635 nm.

Application

Notes:

Additional reported applications (for the relevant formats) include:

immunohistochemistry and fluorescence

microscopy.

Application References:

1. Dolbeare F, et al. 1983. Proc. Natl. Acad. Sci. USA 80:5573.

2. Hirota K, et al. 2007. J. Exp. Med. 204:41.

3. Godebu E, *et al.* 2008. *J. Immunol.* 181:1798.

4. Waskow C, et al. 2008. Nat. Immunol. 9:676.

Description: BrdU is a uridine derivative and a structural analog of thymidine, which can be

incorporated into DNA during the S-phase of a cell cycle as a substitute for thymidine. Cells can be pulse-labeled with BrdU and analyzed with antibodies

BrdU loaded cells, No load control cells

104

10

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products.

Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com

against BrdU to determine the proportion of cells in the S-phase of the cell cycle during a given interval.

Antigen References:

- 1. Barker JM, et al. 2013. PLoS One 8:e63692.
- Duque A and Rakic P. 2011. J. Neurosci. 31:15205.
 Robbins S, et al. 2011. J. Vis. Exp. 55:2855.
- 4. Broekhuizen CA, et al. 2010.