
Product Data Sheet

Alexa Fluor[®] 488 anti-mouse IgD

Catalog # / Size:	2628590 / 100 μg 2628585 / 25 μg
Clone:	11-26c.2a
Isotype:	Rat IgG2a, к
Reactivity:	Mouse
Preparation:	The antibody was purified by affinity chromatography, and conjugated with Alexa Fluor® 488 under optimal conditions.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Concentration:	0.5

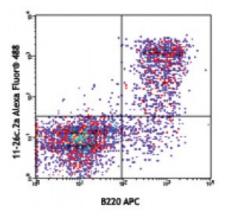

C57BL/6 splenocytes were stained with B220 APC and IgD (11-26c.2a) Alexa Fluor® 488 (top) or rat IgG2a, κ Alexa Fluor® 488 isotype control (bottom).

T" a

Applications:

••		
Applications:	Flow Cytometry	a a
Recommended Usage:	Each lot of this antibody is quality control tested by immunofluorescent staining with flow cytometric analysis. For flow cytometric staining, the suggested use of this reagent is ≤ 0.5 microg per million cells in 100 microL volume. It is recommended that the reagent be titrated for optimal performance for each application.	Lat IgCta Alexa Fluore 488 isotype cont
	* Alexa Fluor® 488 has a maximum emission of 519 nm when it is excited at 488 nm.	
Application Notes:	The 11-26c.2a antibody reacts with immunoglobulin D in all tested mouse haplotypes. The antibody binds membrane IgD expressed on most B cells. The 11-26c.2a antibody neither induces proliferation of splenic B cells nor induces B cell activation. Additional reported applications (for the relevant formats) include: immunohistochemical staining of acetone-fixed frozen sections ^{2,3} .	
Application References:	 Nitschke L, <i>et al.</i> 1993. <i>P. Natl. Acad. Sci. USA</i> 90:1887. (FC) Weih D, <i>et al.</i> 2001. <i>J. Immunol.</i> 167:1909. (IHC) Koni PA, <i>et al.</i> 2001. <i>J. Exp. Med.</i> 193:741. (IHC) Ahuja A, <i>et al.</i> 2007. <i>J. Immunol.</i> 179:3351. (FC) <u>PubMed</u> Haynes NM, <i>et al.</i> 2007. <i>J. Immunol.</i> 179:5099. (FC) Good-Jacobson KL, <i>et al.</i> 2010. <i>Nat. Immunol.</i> 11:535. (FC) <u>PubMed</u> 	

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com 7. Tomayko MM, et al. 2010. J. Immunol. 185:7146. PubMed

- 8. Park SY, et al. 2013. J. Immunol. 190:1094. PubMed
- 9. Rouaud P, et al. 2014. J Exp Med. 211:975. PubMed

Description: Surface IgD is an important B cell differentiation marker.

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com