
APC/Fire[™] 750 anti-mouse lgG1

Catalog # / Size:	2633120 / 100 μg 2633115 / 25 μg
Clone:	RMG1-1
lsotype:	Rat IgG
Immunogen:	Mouse Ig cocktail
Reactivity:	Mouse
Preparation:	The antibody was purified by affinity chromatography and conjugated with APC/Fire™ 750 under optimal conditions.
Formulation:	Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide.
Concentration :	0.2 mg/ml

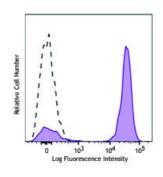

Human peripheral blood lymphocytes were stained with purified anti-human CD3 (clone UCHT1 isotype mouse IgG1, κ) followed by APC/Fire™ 750 antimouse IgG1 antibody (filled histogram) or APC/Fire™ 750 antimouse IgG1 antibody alone (open histogram).

Applications:

Applications:	Flow Cytometry	
Recommended Usage:	Each lot of this antibody is quality control tested. For flow cytometric staining, the suggested use of this reagent is $\leq 0.125 \ \mu g$ per million cells in 100 μ l volume. It is recommended that the reagent be titrated for optimal performance for each application.	
	* APC/Fire™ 750 has a maximum excitation of 650 nm and a maximum emission of 787 nm.	
Application Notes:	S18020H has been tested to block the ligand binding of SCF on c-kit which is determined by the phosphorylation of ERK1/2 in c-kit signaling.	
Application References:	 Rucci F, et al. 2010. P. Natl. Acad. Sci. USA 107:3024. (FC) <u>PubMed</u> Zhang JG, et al. 2011. Int. Immunopharmacol. 11:1685. <u>PubMed</u> 	
Description:	The RMG1-1 monoclonal antibody reacts with immunoglobulin G1 (IgG1) in all tested mouse haplotype (Igh-a and b). It does not react with other isotypes. The RMG1-1 monoclonal antibody may be used as a primary or secondary reagent for ELISA or immunofluorescent analysis.	

For research use only. Not for diagnostic use. Not for resale. Sony Biotechnology Inc. will not be held responsible for patent infringement or other violations that may occur with the use of our products. Sony Biotechnology Inc. 1730 North First Street, San Jose, CA 95112 www.sonybiotechnology.com Antigen References:

- 1. Ledbetter JA, et al. 1979. Immunol. Rev. 47:63. (IHC, IP)
- 2. Hathcock KS. 1991. Current Protocols in Immunology. 3.4.1. (Deplete)
- 3. Takahashi K, et al. 1992. P. Natl. Acad. Sci. USA 89:5557. (Block, IP)
 - 4. Ledbetter JA, et al. 1981. J. Exp. Med. 153:1503. (Block)
 - 5. Hata H, et al. 2004. J. Clin. Invest. 114:582. (IHC)
 - 6. Fan WY, et al. 2001. Exp. Biol. Med. 226:1045. (IHC)
 - 7. Shih FF, et al. 2006. J. Immunol. 176:3438. (FC)
 - 8. Kamimura D, et al. 2006. J. Immunol. 177:306.
 - 9. Bouwer HGA, et al. 2006. P. Natl. Acad. Sci. USA 103:5102. (FC, Deplete)
 - 10. Kao C, et al. 2005. Int. Immunol. 17:1607. PubMed
 - 11. Ko SY, et al. 2005. J. Immunol. 175:3309. (FC) PubMed
 - 12. Rasmussen JW, et al. 2006. Infect. Immun. 74:6590